

Ordering Information

0.014" GW

Balloon Diameter (mm)	Balloon Working Length (mm)				
	15	20	40	80	120
1.5	515-015-41	515-020-41	515-040-41	515-080-41	515-120-41
2.0	520-015-41	520-020-41	520-040-41	520-080-41	520-120-41
2.5	-	525-020-41	525-040-41	525-080-41	525-120-41
3.0	530-015-41	530-020-41	530-040-41	530-080-41	530-120-41
3.5*	-	535-020-41	535-040-41	535-080-41	535-120-41
4.0	540-015-41	540-020-41	540-040-41	540-080-41	540-120-41
4.5*	-	545-020-41	545-040-41	545-080-41	545-120-41
5.0	550-015-41	550-020-41	550-040-41	550-080-41	550-120-41
5.5*	-	555-020-41	555-040-41	555-080-41	-
6.0	560-015-41	560-020-41	560-040-41	560-080-41	-


* available in selected countries, please contact our local sales representative

Compliance Chart

Pressure atm (kPa)	Balloon Diameter (mm)									
	1.5	2.0	2.5	3.0	3.5	4.0	4.5	5.0	5.5	6.0
2 (203)	1.41	1.88	2.37	2.82	3.27	3.70	4.22	4.67	5.09	5.60
4 (405)	1.43	1.90	2.40	2.86	3.32	3.76	4.27	4.74	5.17	5.68
6 (608)	1.45	1.93	2.42	2.89	3.36	3.82	4.33	4.80	5.26	5.76
8 (811)	1.47	1.95	2.45	2.93	3.41	3.88	4.39	4.87	5.34	5.84
10 (1013)	1.48	1.98	2.47	2.96	3.45	3.94	4.44	4.93	5.42	5.92
NOM* 12 (1216)	1.50	2.00	2.50	3.00	3.50	4.00	4.50	5.00	5.50	6.00
14 (1419)	1.52	2.02	2.53	3.04	3.55	4.06	4.56	5.07	5.58	6.08
16 (1621)	1.53	2.05	2.55	3.07	3.59	4.12	4.61	5.13	5.66	6.16
18 (1824)	1.55	2.07	2.58	3.11	3.64	4.18	4.67	5.20	5.74	6.24
20 (2027)	1.57	2.10	2.60	3.14	3.68	4.24	4.73	5.26	5.83	6.32
RBP** 22 (2229)	1.59	2.12	2.63	3.18	3.73	4.30	4.78	5.33	5.91	6.40
24 (2432)	1.60	2.14	2.65	3.22	3.78	4.36	4.84	5.40	-	-
26 (2634)	1.62	2.17	2.68	3.25	3.82	4.42	-	-	-	-

*Nominal Pressure. The nominal in vitro device specifications do not take into account any lesion resistance.

**Rated Burst Pressure - Do Not Exceed (1.5-4.0mm: 22atm, 4.5-5.0mm: 20atm, 5.5-6.0mm: 18atm)

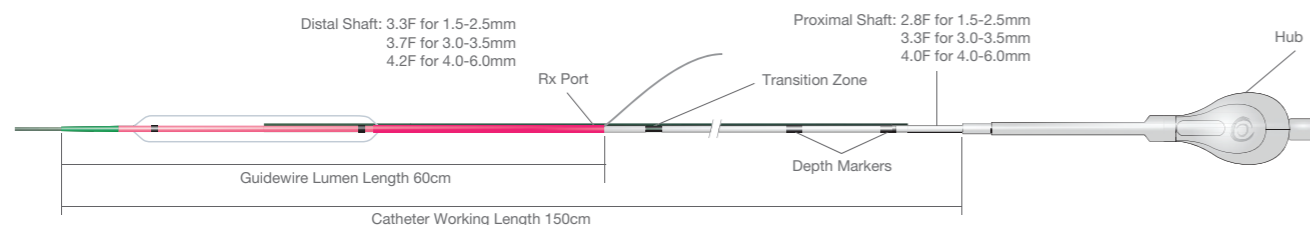

For more information please visit our website at www.OrbusNeich.com or contact us:

EMEA HEADQUARTERS

Drs. W. van Royenstraat 5
3871 AN Hoevelaken
The Netherlands

Phone +31.33.254.1150
Fax +31.33.254.1151

Toll Free Phone 00800.0254.1150*
Toll Free Fax 00800.0254.1151*

* Only for Belgium, Denmark, France, Germany, Ireland, Netherlands, Norway, Sweden and UK

© 2016 OrbusNeich. All rights reserved. **Not available for sale in the USA.**


#G-70-0547 Rev01

JADE and OrbusNeich are registered Trademarks of OrbusNeich Medical, Inc. All data and photos on file. Illustrations are not to be considered as engineering drawings or photos.

For Percutaneous
Transluminal Angioplasty



Strength and precision without compromising deliverability

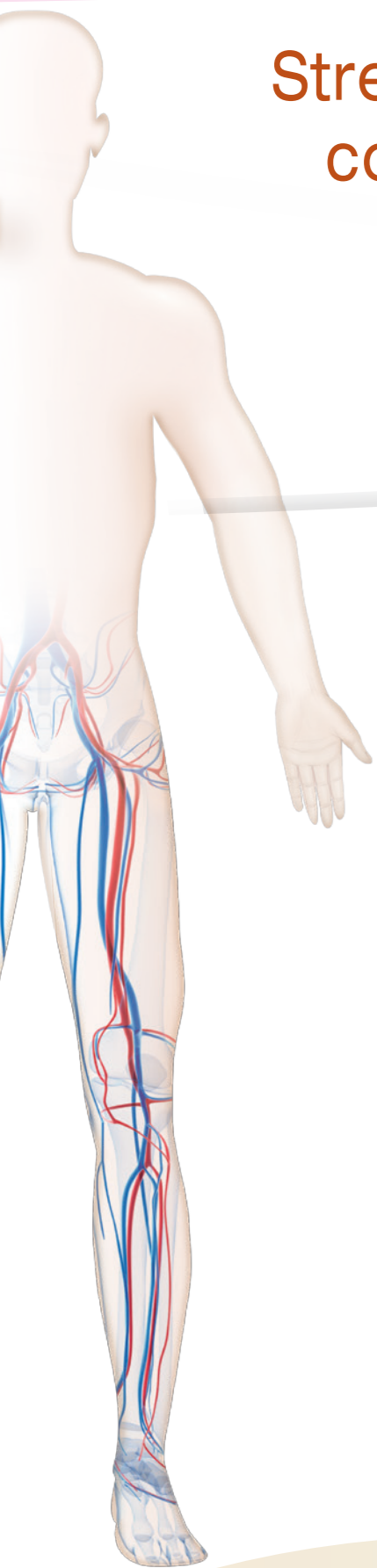

Z-tip technology for managing challenging lesions

High burst pressure **22 RBP** for opening resistant lesions


Hybrid RX/OTW catheter design gives extra pushability for very distal lesions

Technical Specifications

Proximal Shaft	2.8F (1.5-2.5mm); 3.3F (3.0-3.5mm); 4.0F (4.0-6.0mm)
Distal Shaft	3.3F (1.5-2.5mm); 3.7F (3.0-3.5mm); 4.2F (4.0-6.0mm)
Catheter Type	Hybrid RX/OTW
Catheter Working Length	150cm
Tip Length	3.0mm (1.5-3.0mm); 4.0mm (3.5-4.5mm); 5.0mm (5.0-6.0mm)
Guidewire System (GW)	0.014"
Marker Bands	2 marker bands
Nominal Pressure	12atm
Rated Burst Pressure	22atm (1.5-4.0mm); 20atm (4.5-5.0mm); 18atm (5.5-6.0mm)
Coating	Hydrophilic (distal shaft); Silicone (guidewire lumen)
Minimum Sheath Size	4F (1.5-3.5mm); 5F (4.0-6.0mm)

True Non-compliance

Competitor Compliances Comparison for BTK* (2.0mm)
Balloon Diameter (mm)

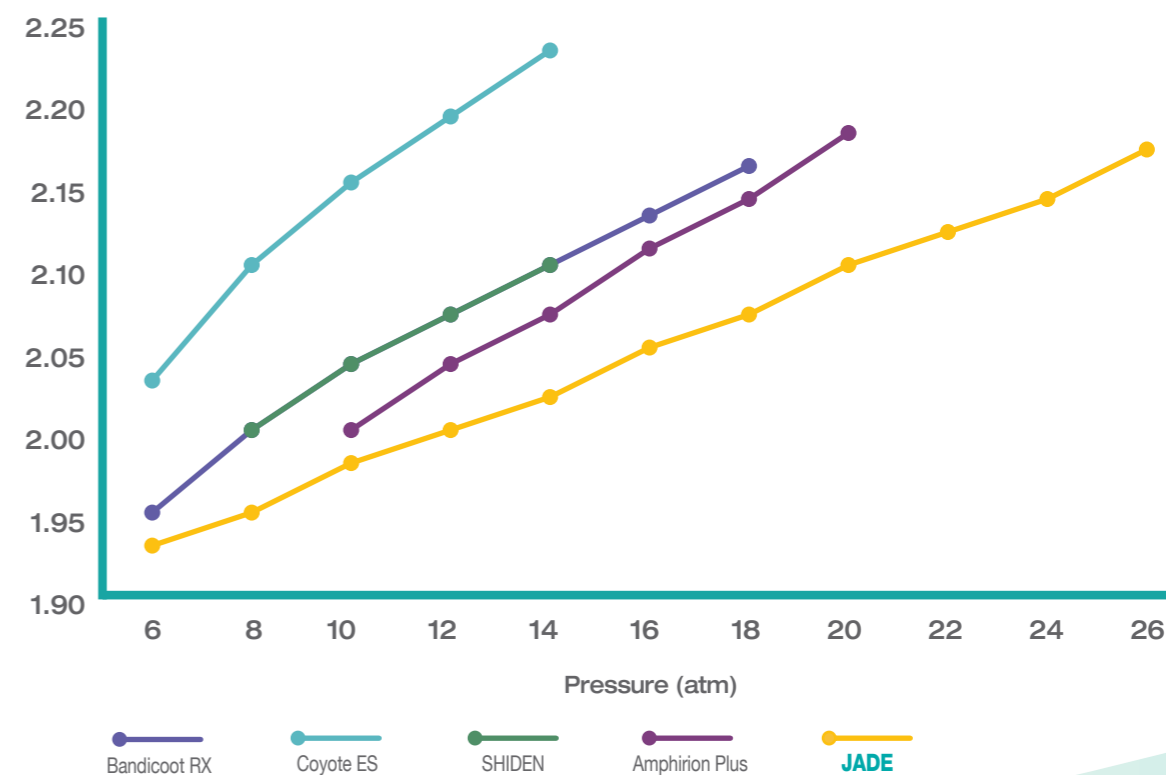

* Based on respective competitor balloons growth rate from NOM to RBP